


Multi-Purpose KNEADING MASSAGER

INSTRUCTION MANUAL

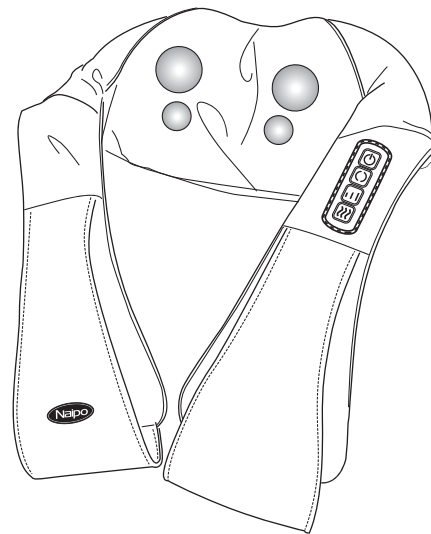

Thank you for purchasing a Naipo Kneading Massager

**PLEASE READ THROUGH THE INSTRUCTIONS
CAREFULLY BEFORE USING THE PRODUCT.**

Shenzhen Bestvap Technology Co., Ltd.

E-mail: bestvap@naipocare.com


Adresa: 503 ManTao Clădirea 720 Minzhi Avenue, Strada Minzhi
District Longhua, ShenZhen, GuangDong 518000



KNEADING MASSAGER

PRODUCT: KNEADING MASSAGER
 MODEL MGS-801
 POWER CONSUMPTION 24W
 VOLTAGE DC 12V 2A
 G.W. 2.1KG

CONTROL PANEL

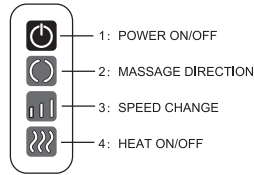

MASSAGER FEATURES

1. Shiatsu-type massage
2. 8 kneading rollers
3. Soothing heat
4. 2 massage directions
5. For neck, shoulders, back, waist, arms, thighs and legs
6. Auto shutdown after 15 minutes
7. Home AC adapter and car adapter included

STORAGE

1. Store the massager in a cool, dry setting.
2. Do not store the massager in the bathroom or in a humid setting with high levels of moisture.
3. Avoid placing the massager near a fireplace, heater, stove, oven or furnace.

TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	SOLUTION
Massager is not turning on	Power plug is loose	Reconnect the plug to the outlet
	Adapter is not plugged in correctly	Unplug the adapter and plug it in again
	Power is off	Push the on button to turn the massager on
	No power	Try a different power outlet or check the car power socket.
Massager is turning off on its own	Power plug is loose	Follow the instructions above and check the connection, adapter, and sockets
off on its own	The 15-minute timer went off	The massager has a timer and will turn off on its own after 15 minutes. You can push the on button to start another 15 minutes session.
Overheat protection is activated during a massage	The massager has an overheat protection feature and will turn itself off when it is overheated.	Please wait for the massager to cool down before powering it back on.

NOTE: IT IS NOT RECOMMENDED TO USE THE MASSAGER FOR MORE THAN 30 CONSECUTIVE MINUTES.

SAFETY INSTRUCTIONS

1. The massager has a built-in timer and it will automatically turn off after 15 minutes. You can push the on button to start another 15 minute session.
 Note: To avoid possible overheating of the unit, it is not recommended to use the massager for more than 30 consecutive minutes.
2. The massager has an overheat protection feature and will turn itself off if it is overheated. Please wait for the massager to cool down before pushing the power button again.
3. Avoid using the massager under direct sunlight or in any high temperature setting.
4. Do not use the massager in the bathroom, bathtub or any setting with high levels of moisture.
5. Disconnect the power when the massager is not being used.
6. Always disconnect the power before cleaning the massager
7. Do not clean the massager with detergent or any sort of liquid.
8. Wipe the massager with a dry or damp cloth to clean the massager.
9. Do not use the massager if any part of the massager (including the power cord) is damaged.
10. Do not let children use the massager.
11. Stop using the massager if you are experiencing discomfort or pain.
12. Do not use the massager if your body is wet.
13. Do not use the massager on sensitive skin.
14. Do not use the massager on joints and bony parts of the body.
15. Do not use the massager on your abdominal region if you are pregnant.



UPPER BACK / NECK



MID BACK



LOWER BACK



THIGHS / LEGS



WAIST



SHOULDERS / ARMS

LIMITED WARRANTY

This product has warranty against defects in materials and workmanship for 90 days from date of purchase. QUALIFICATIONS FOR WARRANTY

Proof of purchase must be presented within 90 days from date of purchase. Defect must not be due to accidental damage or misuse of repair made by unqualified engineers. THIS WARRANTY DOES NOT COVER ANY CLAIMS FOR POTENTIAL DAMAGE OR LOSS. THESE STATEMENTS DO NOT AFFECT YOUR STATUTORY RIGHTS AS A CUSTOMER.