
## User Manual Manual de utilizare

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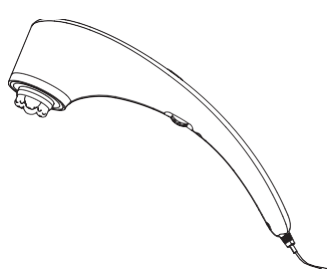#### IMPORTANT SAFETY INSTRUCTIONS

Read all instructions before using the product.  
When using an electrical product, the following basic precautions should always be obeyed:

- ##### DANGER
- To reduce the risk of electric shock:
  - 1) Please disconnect the product from the power supply immediately after usage or before cleaning.
  - 2) Do not reach for the product if it has fallen into water and disconnect it from the power supply immediately.
  - 3) Do not use while bathing or in a shower.
  - 4) Please keep the product away from containers filled with liquids: Please do not use the product in a humid place (like a bathroom).

- ##### WARNING
- To reduce the risk of burns, fire, electric shock, or injuries to people:
  - 1) The product is designed for daily relaxation, not for medical use, or you can use before consulting a doctor.
  - 2) The product should never be left unattended when connected to the power supply. While the product is not in use, or if the parts are to be installed or detached, please disconnect the product from the power supply.
  - 3) Do not use the product under a blanket or a pillow, for excessive heat might occur, and there would exist risks of fire, electric shock, or injury to people.
  - 4) Children or incapacitated people are not recommended to use the product, and close supervision is necessary if they are the users.
  - 5) Use the product only for its intended use as described in this user manual. Use only with attachments recommended by the manufacturer.
  - 6) Do not use the product if the power cable or the plug is damaged. If the product fails to work normally, or if it has dropped into water, please contact the customer service for examination and repair.
  - 7) Do not carry the product by the power cable or use the power cable as a handle.
  - 8) Please keep the product away from hot surfaces.
  - 9) Please do not insert anything into the openings and do not block the air openings during the use to ensure that the product would work normally.

- 01 -

## Contents

## Cuprins

English	01-04
Română	05-08
Français	09-12
Español	13-16
Italiano	17-20

#### Specifications

Dimensions	425*93*130 mm / 16.73*3.66*5.12 inches
Weight	≈1kg
Rated output	15W
Rated voltage	120V- 60Hz

#### Product Operation Guide

- Remove the package.
- Make sure that the product, the power supply, and the power cable are in good condition.
- The product must be turned off when being connected to the product.
- Place the power cable properly and safely to avoid the risk of anyone tripping over it.

#### Massage Areas

#### Q&A

Problem	Cause	Solution
The product does not work.	The product is not connected to the power supply. The product has not been turned on.	Please connect the product to the power supply and turn on it. Please turn on the product.
The product works noisily.	The massage node has not been installed properly.	Please install the massage node again.

- 03 -

#### Product Maintenance

- Disconnect the product from the power supply before cleaning it.
- Only clean the product in the specified manner. Avoid any liquid entering inside the product or accessories.
- Do not use the product again until it is completely dry.

#### Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us via the email or the address below.  
Web: www.naipocare.com  
E-mail: support@naipocare.com  
Address: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China.

\*Note: Products purchased directly from Naipo benefit from our after-sales service, if the purchases are made from other selling terms, we recommend you contact them for warranty issues.

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