FIND WELL-BEING

# USER MANUAL

## CONTACT US

Aukey International Limited

Web: [www.naipocare.com](http://www.naipocare.com)

Email: [wholesale@naipocare.com](mailto:wholesale@naipocare.com)

Add: Room 601 6/F, Tung Hip Commercial Building, 248 DesVoeux  
Road Central, Sheung Wan, Hong Kong

## MADE IN CHINA

[/naipocare](#) [/naipocare](#) [/naipocare](#)

**oCuddle™ Plus**

SHOULDER MASSAGER

# CONTENTS |

English

02

---

Româna

10

---

Français

20

---

Español

30

---

Italiano

40

---

## IMPORTANT SAFETY INSTRUCTIONS

When using the product, all the following instructions and basic precautions should be taken.

## DANGER

To reduce the risk of electric shock, burns, fire, etc.:  
Always unplug this product from the electrical outlet immediately after use and before cleaning.  
Please do not disassemble the product by yourself.

## WARNING

### Danger - To reduce risk of electric shock:

Do not operate under a blanket or a pillow and allow normal heat dissipation, otherwise the massager could overheat, which might cause fire, electric shock, or injury.

Do not unplug by pulling the power cord to prevent damage.

Only use this product with the power supply unit provided, and only with the voltage specified on the label.

Do not operate in places where corrosive products are used or where oxygen is administered.

Do not use in environments with high humidity (eg: bathroom) or high temperature, which could lead to electric shock or other safety hazards.

There will be high temperature in the car, please do not put the product in the car.

Do not use outdoors.

### Attention - Read all instructions carefully before using the product:

This product is a personal massage product, not a medical instrument, and it cannot be used for other commercial activities.

Only use this product for the intended purpose as described in this manual.

Please switch the power off when the product is not in use.

The product should never be left unattended when plugged in. Unplug from outlet when not in use.

The product is equipped with a heating function. People insensitive to heat should be careful when using.

Keep the cord away from heated surfaces, so as not to cause serious burns.

Please pay special attention to high wear components such as the massage heads. Replace defective components immediately and/or keep the product out of use until repair.

Do not use the product if it has a damaged cord or plug, if it is not working normally, if it has been dropped or damaged, or if it has dropped into water.

Stop using the product immediately if you feel discomfort or pain during use.

Stop using the product if it is not working normally.

Do not use the product while driving, and observe the safety driving regulations.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision. Cleaning and user maintenance shall not be made by children without supervision.

**If you have any of the following conditions, please consult doctor before using this product:**

People with physical disabilities.

Those with chronic illnesses.

Those who are pregnant.

Those who are menstruating.


Those with physiological issues.

People with pacemakers or concerns about their health.

## MASSAGE AREAS
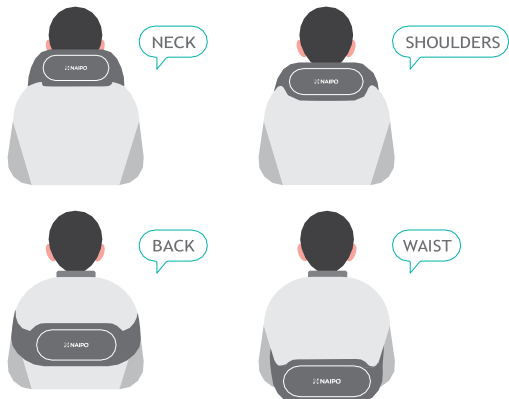
## IN THE BOX

Massager


Wall Adapter

## OPERATION GUIDE

### How to use

- Place the product on desired place.
- Long press (2 seconds) the power button to start the product. The heating and rollers will automatically switch on. The default timer is 15 minutes.
- Press the heating button to switch between 2 levels (High / Low) heating and ON / OFF.
- Press the massage direction adjustment button to change the massage direction.
- Press the speed adjustment button to choose among 3 massage speed levels (Low / Medium / High).
- Press the power button to stop the massage at any time.
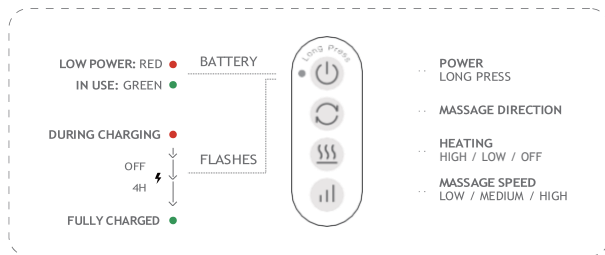
### Charging

- The indicator will light up in red when the product is low on battery. When the battery is low, plug in and charge.
- Insert the cylindrical end of the charger cable into the product charging port, and plug the adaptor into an electrical power outlet.
- A full charge cycle takes approximately 4 hours. During charging, the indicator flashes red and it turns to green when the product is fully charged.
- When connected to a power supply, the product will automatically charge if the battery is low.
- The product can be used continuously for 60 minutes in a fully-charged condition.

**NOTE:** The product can be used while plugged into a power source but the product will not be charged.

## CLEANING AND STORAGE

To prevent damage and discoloration, do not clean the product with corrosive detergent. Do not write on, or poke the product with sharp pen tips.

Before cleaning, make sure that the product has been turned off and that the adapter has been pulled out from the outlet. Clean with a dry cloth.

The zippered fabric cover is removable for cleaning.

Keep the product clean and store in a safe place. Keep away from environments with high humidity or temperature.

## MALFUNCTIONING AND MAINTENANCE

- If the product cannot be started:
  - (1) Fully charge the product, which takes 4 hours.
- If the product is unable to charge:
  - (1) Reconnect the charger to the power port;
  - (2) Try a different power socket.
- If the product stops during use:
  - (1) Auto timer has gone off;
  - (2) It has stopped due to excess pressure. Let it rest before restarting.
- If the product appears damaged:
  - (1) Do not continue to use. Consult dealer or technical personnel with maintenance experience for help.

**NOTE:** Do not disassemble or modify the product by yourself. NAIPO will not be responsible for any loss caused by unauthorised disassembly or modification of the product. This action may put the user in danger and/or render the product unusable. Please contact NAIPO customer service for further assistance.

**NAIPO shall not be liable for costs, damages or repairs incurred as a result of, due to, or from:**

- The warranty period having expired.
- Unauthorized actions such as: product repair, misuse, abuse, modification, disassembly, etc.
- Artificial or accidental damage.
- Damaged internal and external packaging or missing accessories.
- Damage caused by force majeure.

**This product has an internal, non-replaceable battery. This battery is not user-serviceable. Please dispose in accordance with local, state, province, and country regulations.**

## SPECIFICATIONS

PRODUCT NAME	oCuddle™ Plus Shoulder Massager
MODEL	oCuddle- P1
NET WEIGHT	1.8kg   3.9lb
DIMENSION	430*150*200mm   17.0*6.0*8.0in
BATTERY	2200mAh
INPUT	14V --- 1.7A
RATED POWER	23.8W
BRITISH STANDARD ADAPTER NUMBER	Use only with power supply unit AP049B-14170
EUROPEAN STANDARD ADAPTER NUMBER	Use only with power supply unit AP049E-14170

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